

## Anti-Dis2 p-T316 (S. pombe) antibody, rabbit serum

63-121 50  $\mu$  l

**Key words:** Phoshoprotein phosphatase, DNA damage checkpoint, microtubule cytoskeleton organization, rRNA processing, signal transduction, signal transduction, mitotic cell cycle, homologous chromosome segregation

**Backgroud:** S. pombe dis2 gene encodes Serine/threonine-protein phosphatase PP1-1 (327 aa, 37.6 kDa) which plays essential role in cell cycle control and required for exit from mitosis.

Dis2 protein is phosphorylated at Thr 316 in mitosis

## **Applications**

1. Western blotting (1/1,000~1/2,000 dilution). Not tested for other applications

Immunogen: Synthetic peptide NWHMT(PO3)PPRKN conjugated to KLF

Specificity: Reacts with S. pombe Dis2 protein phosphorylated at Thr316. Not tested with other species

Form: Rabbit antiserum added with 0.05 % sodium azide

Storage: Shipped at 4°C and upon arrival, centrifuge briefly, aliquot and store at -20°C.


Database Links: <u>UniProt/Swiss-Prot P13681</u>

## PomBase SPBC776.02c

**References:** This antibody was described in Ref. 1 and used in the following publications.

- 1. Ishii K. et al (1996) Requirement for PP1 phosphatase and 20S cyclosome/APC for the onset of anaphase is lessened by the dosage increase of a novel gene sds23+. EMBO J. 15: 6629-40. PubMed 8978689 WB
- 2. Sutani T. et al (1999) Fission yeast condensin complex: essential roles of non-SMC subunits for condensation and Cdc2 phosphorylation of Cut3/SMC4. Genes Dev. 13: 2271-83. PubMed 10485849

  WB

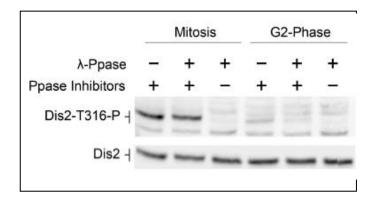

Related Product: 63-119 anti-Dis2 antibody

Identification of Dis2 Figure. phosphorylated at T316 by western blotting with the antibody. S. pombe crude extracts prepared from mitotic and G2-phase were analysed by WB. Dis2 and phosphorylated Dis2 were detected with anti-Dis2 antibody (BA 63-119) and this antibody, respectively. Phosphorylation is increased in mitosis and sensitive to  $\lambda$ -phosphatase. The antibody was used at 1/1,000 dilution in PBS containing 0.1% Tween and 1% milk. Courtesy of Mr M. Swaffer at Cancer Research UK

BioAcademia, Inc. Tel. 81-6-6877-2335 Fax. 81-6-6877-2336 info@bioacademia.co.jp http://www.bioacademia.co.jp/en/